
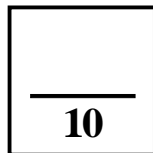


Science 10-Biology

Activity 5

Observing Plant and Animal Cells

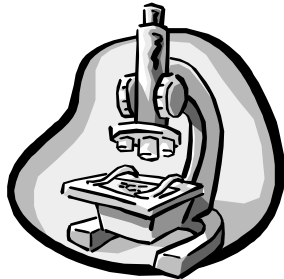

Name \_\_\_\_\_

Due Date \_\_\_\_\_

Show Me  Hand In

**Correct and Hand In Again By** \_\_\_\_\_

**Purpose:**

- To observe a prepared slide of plant cells and a prepared slide of animal cells.
- To locate and draw some of the organelles in these cells.
- To prepare a wet mount slide of living plant cells, stain it and observe it.

**Part 1-Prepared Plant Cell Slide**

**Procedure:**

1. Take a prepared slide of plant cells (spirogyra, leaf cross section, etc.) . First focus it under Low power, center it then focus under Medium power. Center it again and focus as well as you can under High power. (Magnification = 400X)  
As you look at the image, answer the following questions:

Estimate how many cells in total you see in your field of view? \_\_\_\_\_ cells

Estimate how many different kinds of cells you see. \_\_\_\_\_ kinds


On the diagram *on the next page*, **sketch** several cells (an example of each type you see) and **label** the following organelles:

- Cell Wall**
- Cell Membrane**
- Central Vacuole**
- Nucleus**
- Chloroplast**

Use the diagram on **page 332 of the Text** to help you identify the parts you are looking at.



Slide Name: \_\_\_\_\_

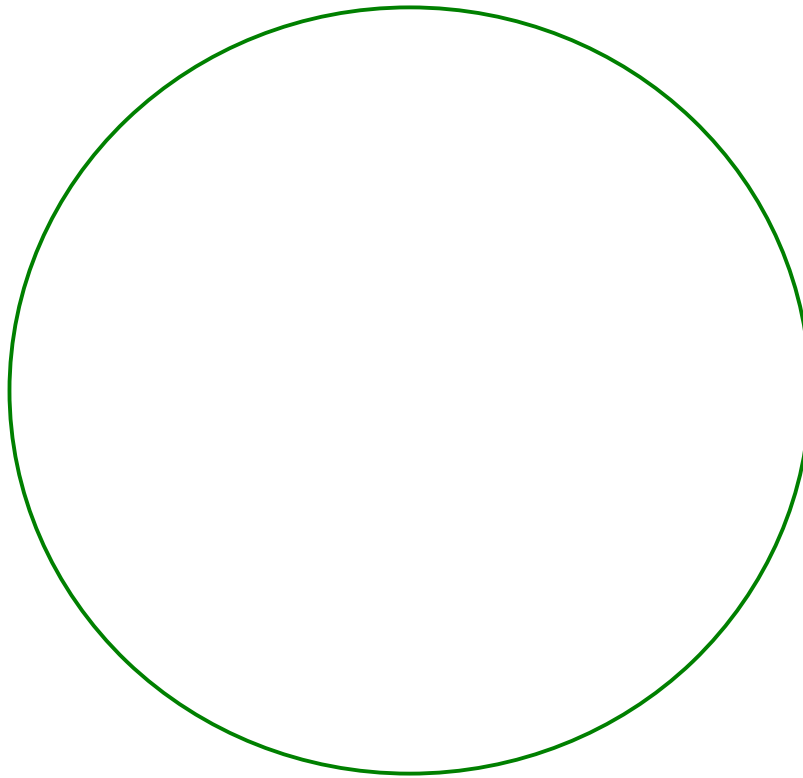

Microscope Magnification \_\_\_\_\_ X

2. State the **function** of each of the following:

**Cell Wall** –

**Chloroplast** –

**Vacuole** –

**Nucleus** –

**Cell Membrane** –

Why do you think it is hard to see the **cell membrane** in plant cells?

## Part 2-Prepared Animal Cell Slide


### Procedure:

1. Take a prepared slide of an animal cell (epithelium, liver, etc.)  
Animal cells tend to be very specialized and it is not easy to get a “typical” animal cell. It is also not easy to distinguish animal cells that are side by side.

Why do you think this is often the case?

On the diagram below, **sketch** several cells and **label** the following parts:  
**Cell Membrane, Nucleus, Cytoplasm**

Slide Name: \_\_\_\_\_

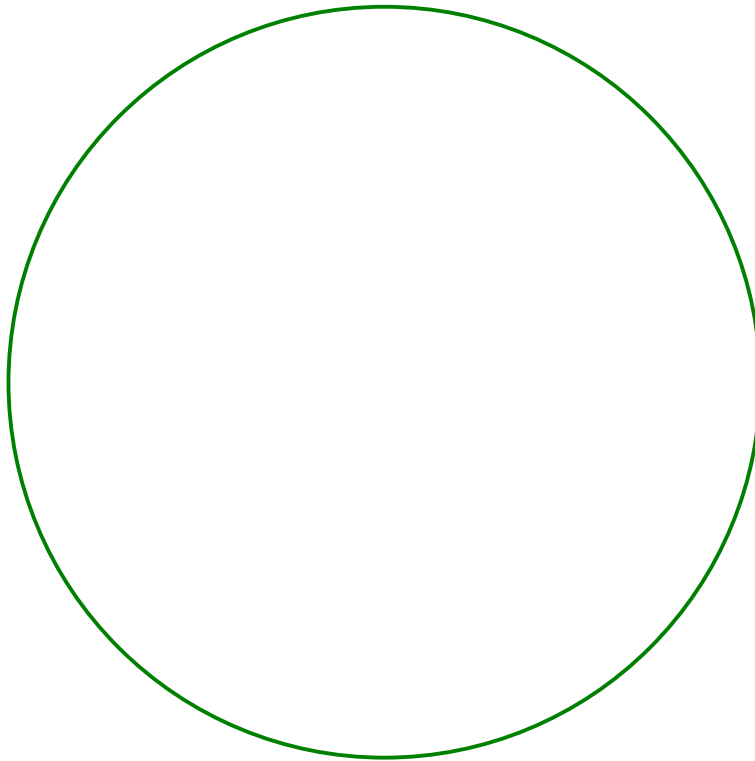


Microscope Magnification \_\_\_\_\_X

### Part 3-Living Plant Cells

#### Procedure:

- Using a scalpel, scissors, tweezers etc. pull a **thin** film of tissue from a piece of onion. It must be very thin and fairly transparent. Make a wet mount (See the diagram on the left of page 320 in the Text.)  
Focus on Low power first, get a good view, then move up to Medium or High power. Sketch what you see and label any parts you can recognize (cell walls, cell membranes, vacuoles, nuclei, chloroplasts etc.)

Slide Name: \_\_\_\_\_

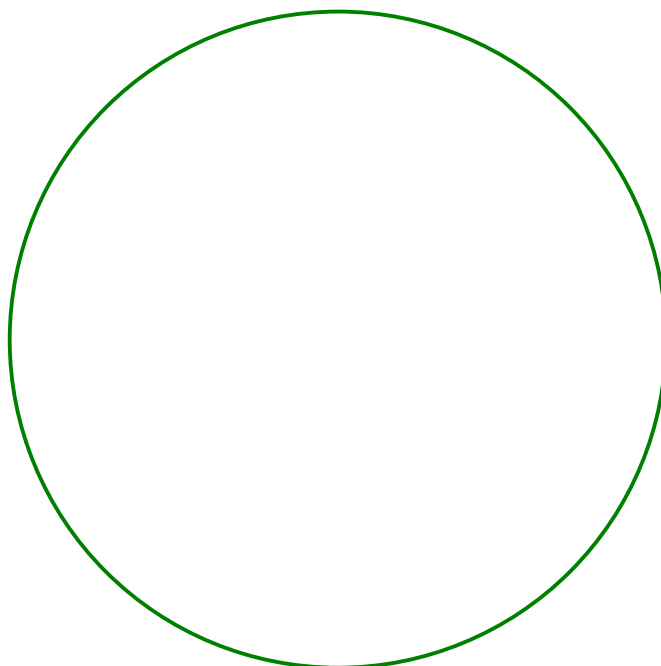

Microscope Magnification \_\_\_\_\_X

- Remove the slide from the microscope. Using a piece of paper towel, draw out the water from under the cover slip. (See the diagram on the left of page 333 in the Text). Place a drop or two of **IKI solution** so that it goes **under the cover slip**. Put your slide back under the microscope and focus on Medium or High power again. ***IKI is a stain that turns starch a blackish colour.***

Do you see any evidence that there is **starch** in onion cells? \_\_\_\_\_

What parts **stained** (went blackish)? \_\_\_\_\_

What parts **did not stain**? \_\_\_\_\_

How does staining cells help us to see them?